

PRODUCT HIGHLIGHTS:

- DURABLE, LIGHTWEIGHT UNIBODY PLATFORM
- LOW BUILD HEIGHT/HIGH DISPLACEMENT
- RESPONSIVE SHOCK ABSORPTION
- INCREASED GROUND COMPLIANCE
- FLUIDITY THROUGH THE GAIT CYCLE
- 40° TORSIONAL ROTATION
- CORROSION RESISTANT AND WATERPROOF

OVERVIEW:

The Fiberglass META Shock X delivers the full benefits of the META-Unibody® platform in fiberglass for K3/K4 prosthetic users up to 365 lbs. The split keel design provides users improved stability and ground compliance, enabling conformance and easy adaptation to variations in terrain, while offering a smooth rollover experience.

The revolutionary unibody design and low build height, characteristic of the META foot family, offer market-leading levels of deflection in a compact, lightweight package. The Fiberglass META Shock X builds upon this foundation by combining multiphase shock absorption with 40° axial rotation in the ankle.

The integrated torsion component promotes natural rotational movement and gait fluency with less torque applied to the residual limb, delivering patient comfort via reduced stress on the residual limb.



PDAC
APPROVED



✓ TRANSFEMORAL ✓ TRANSTIBIAL

PDAC APPROVED L-CODES: L5987 and L5984

FIBERGLASS META SHOCK X: SPECIFICATIONS


HEEL HEIGHT: 10 mm (3/8")

BUILD HEIGHT: 138 mm (5.43")

**APPROXIMATE
PRODUCT WEIGHT:** 815 g

*Std. size 26 Cat. 5 with foot shell and sock

TRIAL PERIOD: 60 days
from date of patient fitting

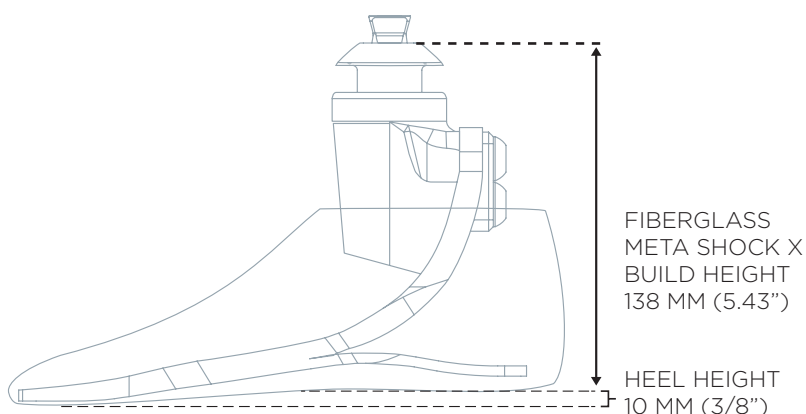

FIBERGLASS META SHOCK X: ORDERING INFORMATION

FIBERGLASS META SHOCK X CATEGORY CHART				
Adjusted User Weight (patient body weight + loads carried)	Standard (Low to Medium) Impact Level		High Impact Level	
	Sizes 22-24 cm	Sizes 25-30 cm	Sizes 22-24 cm	Sizes 25-30 cm
99-115 lbs (45-51 kg)	1	1	1	1
116-130 lbs (52-58 kg)	1	1	2	2
131-150 lbs (59-67 kg)	2	2	3	3
151-170 lbs (68-76 kg)	3	3	4	4
171-194 lbs (77-87 kg)	4	4	5	5
195-221 lbs (88-99 kg)	5	5	-	6
222-256 lbs (100-115 kg)	-	6	-	7
257-287 lbs (116-130 kg)	-	7	-	8
288-324 lbs (131-147 kg)	-	8	-	9
325-365 lbs (148-165 kg)	-	9	-	-

**Note: The Standard impact level is recommended for most patients for day-to-day use; high impact level should be chosen for active outdoor users engaged in activities that involve heavy ground impact forces (e.g., sports, construction work)*

HOW TO ORDER							
Fiberglass META Shock X Foot with Spectra Sock				Sandal Toe Foot Shell			
Ordering example: Fiberglass META Shock X Foot, 25 cm, left foot, toe resistance 2, medium impact				Ordering example: Sandal Toe Foot Shell, 25 cm, left foot, tan			
M006 - 25L - 2M				FS01 - 25LT			
Model	Size of Foot	Foot	Toe Resistance	Impact Level	Size of Foot	Foot	Color
M006 = Fiberglass META Shock X Foot	22-30 cm	L = left R = right	01-09 (Refer to the chart above)	L = low M = medium H = high (See the note above)	22-30 cm	L = left R = right	B = Buff T = Tan M = Medium Brown

FOR MORE INFORMATION CALL 800.848.4930 OR 740.869.3377
OR VISIT WILLOWWOOD.COM

WillowWood®